



BEFORE & AFTER TREATMENT




BEFORE TREATMENT

- 15mm VP on 1st MPT Joint
- 2 treatments of 5W for 5 seconds
- Resolved within 12 weeks
- 6 month follow-up showed no scarring & no recurrence*



AFTER 3 WEEKS

(3 weeks after 1st treatment)

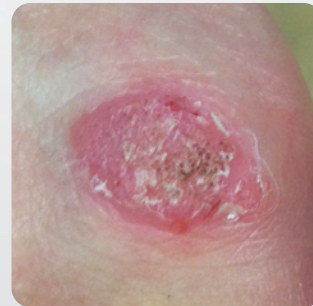

AFTER 9 WEEKS

(6 weeks after 2nd treatment;
no further treatment required)



AFTER 12 WEEKS


(No further treatment after
2nd treatment)

For more information visit www.treatwithswift.com
or speak to your Podiatrist today




MICROWAVE THERAPY FOR SKIN


Microwave energy is widely used in a number of medical fields to elevate tissue temperatures and create precise, localised cell destruction.

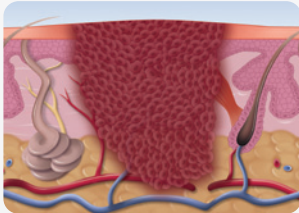

Infected tissue can exist several millimetres below the surface and can often be difficult to treat using traditional methods, resulting in either untreated tissue or significant damage.

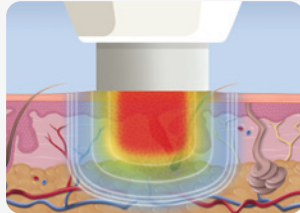

Swift delivers a precise, highly controlled energy dose. As microwaves travel into the tissue, water molecules begin colliding and creating localised heat energy – quickly destroying all infected tissue within a predetermined depth.

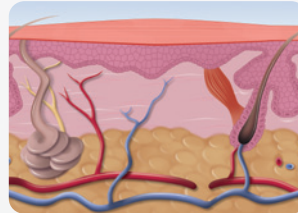

In just seconds the treatment is complete, leaving the body to absorb and replace the treated tissue.

For more information visit www.treatwithswift.com
or speak to your Podiatrist today

